

eGateway

INTEGRATION SOLUTION

Features and Benefits

- Supports vital sign monitors, networked bedside monitors, central stations and A-Series anesthesia systems
- Provides “single pipe” for data interchange between Mindray monitoring solutions and hospital information systems
- Supports HL7 versions 2.1 – 2.6 including IHE PCD profile to speed integration
- Enables download of patient demographics from two different ADT servers via optional ADT interface to simplify admission and improve accuracy, supporting positive patient identification
- Transfers vital signs and physiologic data to as many as four hospital data systems (HIS, CIS, EMR, etc.) with standard Results Reporting interface
- Spot Check interface supports Accutorr 7 vital sign monitors; provides unsolicited communication to EMRs including manual parameters and modifiers identified at bedside
- Alarm information output via alarm interface to as many as three alarm management platforms; sends alarms and waveform snippets using IHE PCD framework
- Automatically detects devices on network, providing up to 48 hour backfill to hospital EMR
- Employs Network Time Protocol (NTP) to synchronize time across Mindray solutions with hospital network time



Providing meaningful interoperability between Mindray clinical solutions and hospital systems


The **Mindray eGateway Integration Solution** is a powerful informatics component for bi-directional communications between networked Mindray clinical solutions and hospital ADT, CIS, EMR and alarm management systems. Interoperability is simplified through the adoption of standard HL7 protocols and the IHE PCD profile for seamless connectivity. For flexibility, the eGateway is available as a hardware and software solution or as a software only solution capable of VMware deployment.

Effectively managing patient demographics, EMR results data, alarm and waveform distribution, the eGateway can support the entire patient monitoring network across the healthcare enterprise.

mindray

eGateway Network Diagram

The eGateway bridges the Mindray Clinical Network and the Hospital Network.

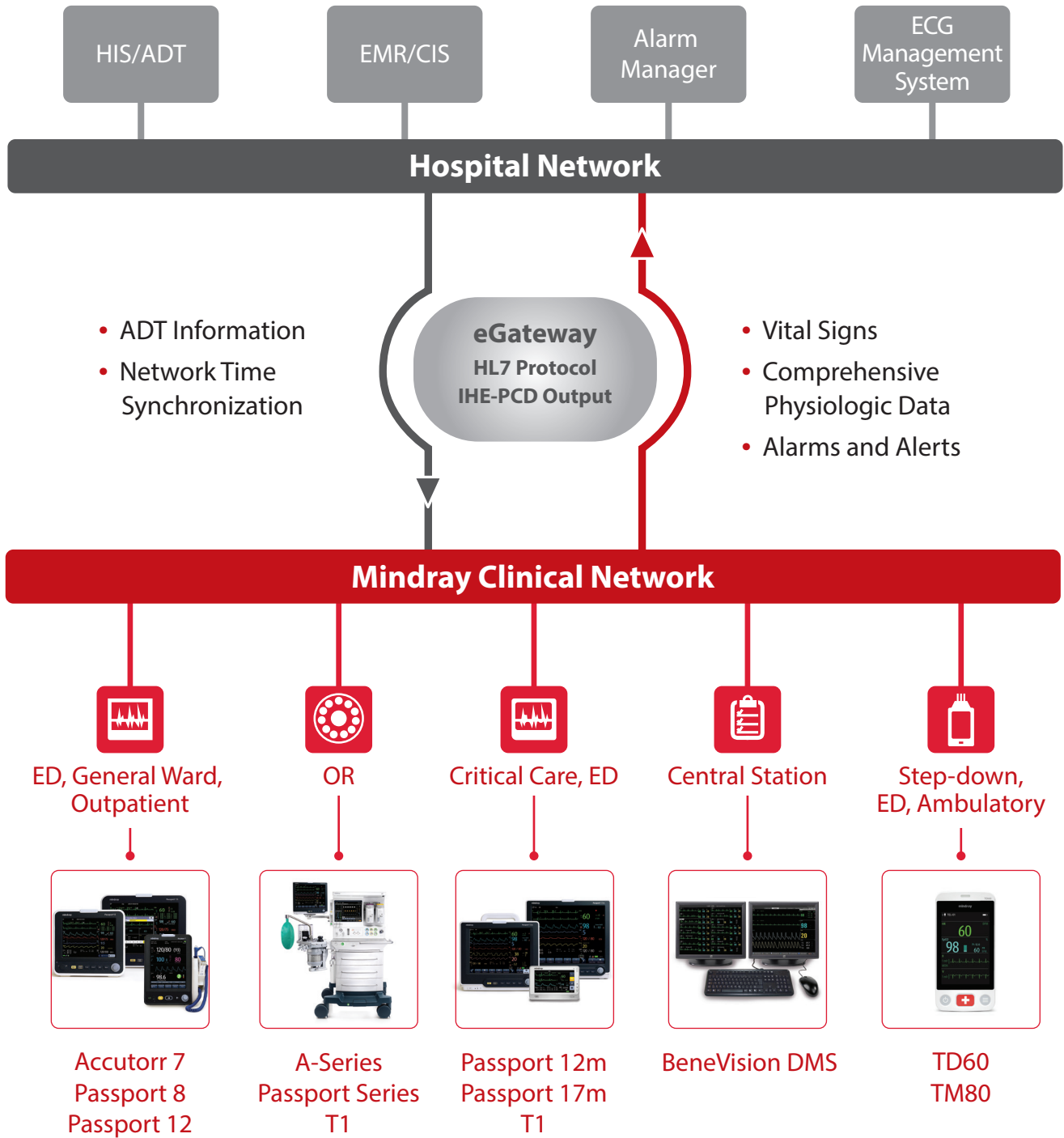

Mindray North America
800 MacArthur Blvd., Mahwah, NJ USA 07430
Tel: 800.288.2121 or 201.995.8000 Fax: 800.926.4275 www.mindray.com

©2018 Mindray DS USA, Inc. Mindray® is a registered trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. All brands and product names are trademarks of their respective owners. Specifications subject to change without notice. Contact your Mindray Representative for the most current information. 07/18 P/N: 0002-08-7009 Rev E

mindray
NORTH AMERICA