

eGateway

INTEGRATION SOLUTION

Features and Benefits

- Supports vital sign monitors, networked bedside monitors, the A-Series anesthesia systems and the BeneVision DMS
- Provides “single pipe” for data interchange between Mindray monitoring solutions and hospital information systems
- Supports HL7 versions 2.1 – 2.6 including IHE PCD profile
- Enables download of patient demographics from ADT servers via optional ADT interface to simplify admission and improve accuracy, supporting positive patient identification
- Document sharing channel for convenient e-file of real-time and historic data from the DMS to the EMR or other destination replacing paper printing
- Spot Check interface supports Accutorr 7 vital sign monitors; provides unsolicited communication to EMRs including manual parameters and modifiers identified at bedside
- Alarm information output via alarm interface to as many as three alarm management platforms; sends alarms and waveform snippets using IHE PCD ACM framework
- Supports receiving up to 48 hours of backfill data from patient connected devices and backfills data to the EMR when necessary
- Employs Network Time Protocol (NTP) to synchronize time across Mindray solutions with hospital network time
- High Resolution Interface for real-time distribution of patient parameters (PCD), waveforms (WCM) and alarms (ACM)



Providing meaningful interoperability between Mindray clinical solutions and the Electronic Health Record


The Mindray eGateway Integration Solution is a powerful informatics component for bi-directional communications between networked Mindray clinical solutions and hospital ADT, CIS, EMR and alarm management systems. Interoperability is simplified through the adoption of standard HL7 protocols and the IHE PCD profile for seamless connectivity. For flexibility, the eGateway is available as a hardware and software solution or as a software only solution capable of VMware deployment.

Effectively managing patient demographics, EMR results data, alarm and waveform distribution, the eGateway can support the entire patient monitoring network across the healthcare enterprise.

mindray

eGateway Network Diagram

The eGateway bridges the Mindray Clinical Network and the Hospital Network.

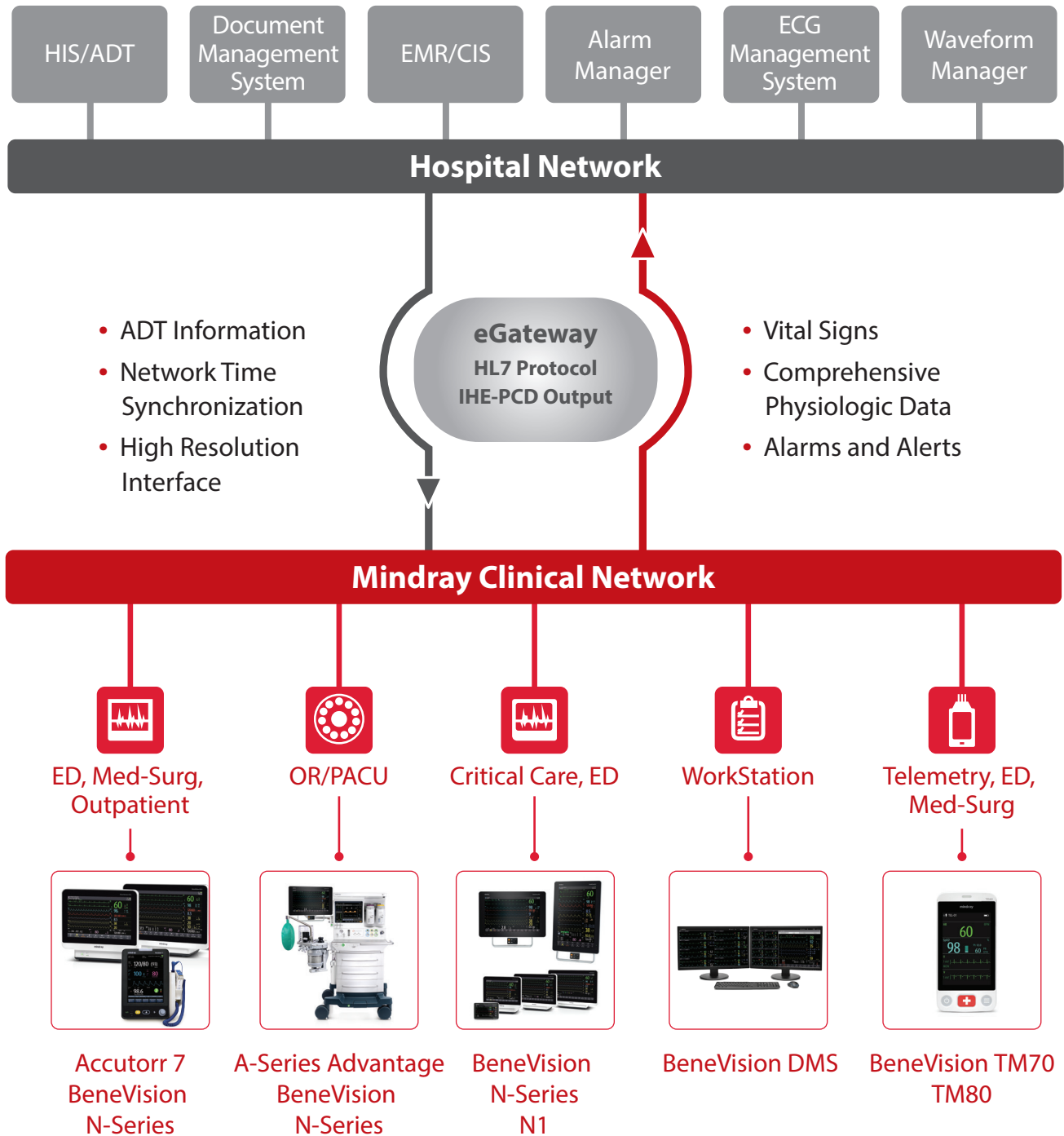

Mindray North America
800 MacArthur Blvd., Mahwah, NJ USA 07430
Tel: 800.288.2121 or 201.995.8000 Fax: 800.926.4275 www.mindray.com

©2021 Mindray DS USA, Inc. Mindray® is a registered trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. All brands and product names are trademarks of their respective owners. Specifications subject to change without notice. Contact your Mindray Representative for the most current information. 01/21 P/N: 0002-08-7009 Rev G

mindray
NORTH AMERICA