

Accutorr[®] 7

VITAL SIGNS MONITOR

Features and Benefits

- Large 8.4" back-lit high resolution TFT display provides clear view
- Intuitive touchscreen offers easy access to information
- Manual, automatic and customizable NIBP measurement modes offer flexibility for various clinical applications
- Available SpO₂ technology includes Masimo SET[®] (standard) or Nellcor[®] Oximax[®]
- PI (Perfusion Index) from SpO₂ can guide caregivers to best measurement location
- Temperature options include: Exergen Temporal Scanner[™] and Mindray SmarTemp[™]
- Low flow sidestream CO₂ provides measurements of FiCO₂, etCO₂, awRR with a CO₂ waveform
- Standard MEWS (Modified Early Warning Score) offers effective support for clinical decisions
- Manual input of vitals such as respiration rate, pain level, glucose, fluids and degree of consciousness
- Stores up to 5000 measurements, accessible by patient ID for easy review
- Choice of hardwired or wireless EMR connectivity options
- ADT and results messaging easily accomplished via Mindray's powerful eGateway integration solution



The Accutorr 7 Vital Signs Monitor

Meeting the needs of low acuity areas, the Accutorr 7 is an easy to use vital signs monitor that is ready to go when you are. From spot check to continuous monitoring, assess patients' vitals with NIBP, SpO₂ and temperature measurements. Exergen Temporal Scanner[™] Thermometer quickly measures patient temperature and provides significant cost savings by eliminating probe cover usage. With the addition of low flow sidestream CO₂ and central monitoring connectivity, identify signs of patient deterioration due to respiratory distress. The integrated early warning scoring system can aid in clinical decisions. With its flexibility, the Accutorr 7 is a vital solution for delivering quality care.

The Accutorr 7 is designed to be connected. In today's healthcare environment, vital sign data collection, validation and timely EMR posting are essential in providing the best possible medical care. The Accutorr 7 offers flexible connectivity options meeting your institutional requirements. With wireless capabilities using shared infrastructure, mobility between patients and real time data posting at point-of-care to your EMR can improve efficiency and reduce manual errors.

mindray

Accutorr® 7

Connect. Capture. Send.


Seamlessly connect Accutorr 7 to your hospital's ADT and EMR system with the capability for clinician authentication to help simplify workflow, reduce errors and manage time.



- At point-of-care, clinicians can validate patient information and automatically send vitals to the patient chart, improving efficiency
- Optional medical grade barcode scanner can speed patient ID input, supporting 1D and 2D barcodes

Designed to be Flexible when Connected

- Wireless and wired Accutorr 7 connectivity options available
- Centrally monitor patients on the Accutorr 7 with the Mindray BeneVision Central Station
- Connect the Accutorr 7 with ADT and EMR systems via Mindray's powerful eGateway Integration Solution
- Clinician authentication through eGateway MLDAP application connecting with the hospital Active Directory (AD)
- Direct HL7 interface capabilities
- Standard RS232 interface for third party connectivity solutions

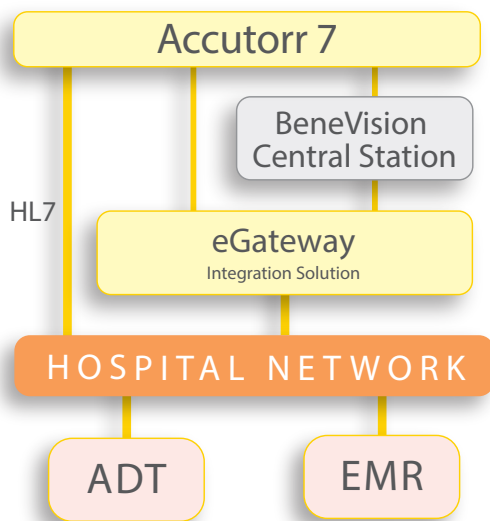


Recognize respiratory distress and quickly respond to signs of patient deterioration.

Mindray's low flow sidestream CO₂ technology with a sampling rate as low as 50 ml/min meets measurement requirements for patients with low tidal volume.

The Accutorr 7 with CO₂ is compatible with Covidien® Microstream™ FilterLines™.

Clinical Decision Support with Modified Early Warning Score (MEWS)

Standard MEWS offers support for clinical decision making with customizable parameters based on facility protocols.

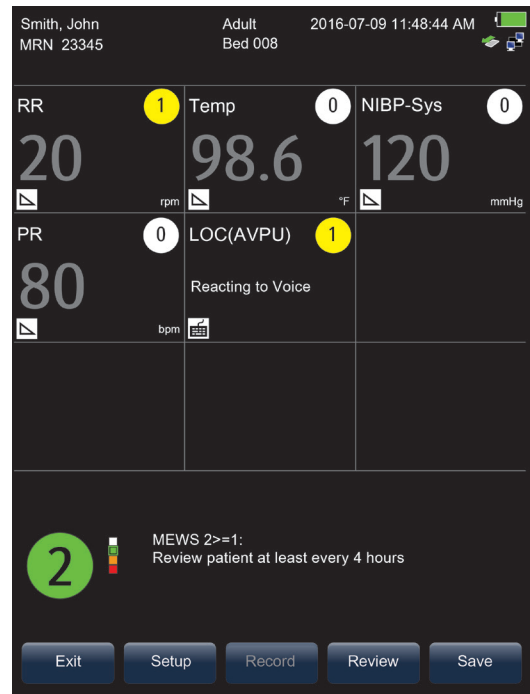

Mindray North America
 800 MacArthur Blvd., Mahwah, NJ 07430
 Tel: 800.288.2121 Tel: 201.995.8000 Fax: 800.926.4275 www.mindray.com

©2020 Mindray DS USA, Inc. Mindray® is a registered trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. All brands and product names are trademarks of their respective owners. Specifications subject to change without notice. Contact your Mindray Representative for the most current information. 04/20 P/N: 0002-08-6373 Rev G

